


# 4D Clean Kit

Quick Start Guide

Version: V1.0.40828


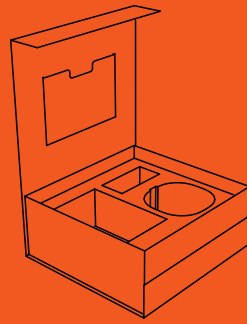


4D Clean Kit rejuvenates your 4D Aligner™ daily to ensure best possible performance and aesthetics.


## 01 Box Contents



1 4D Clean Kit

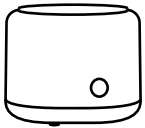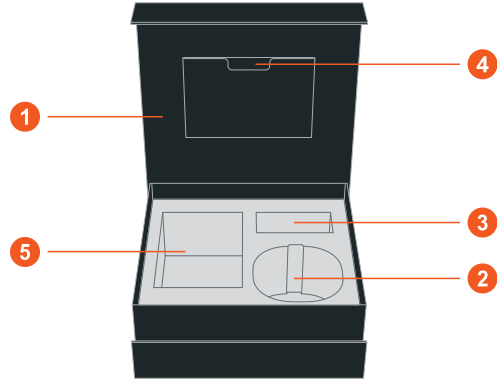

2 Ultrasonic Cleaner x 1

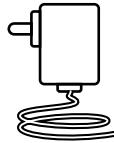

3 Power Adapter x 1




4 Quick Start Guide




5 Cleaning Tablets x 96


## 02 4D Ultrasonic Cleaner

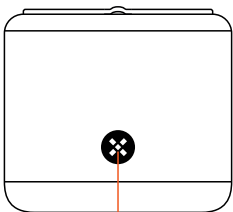

On/Off Button  
Indication Light

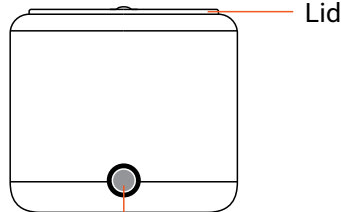

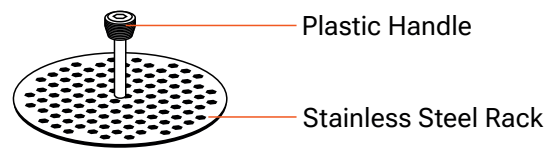

Power Input



Filter Rack


### Light Indication:

 **Black:**  
Adapter not connected

 **White:**  
Adapter connected but power off

 **Orange:**  
Power on

## 03 Operating 4D Ultrasonic Cleaner

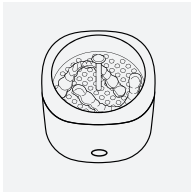


### Step 1:

- Place appliance(s) side by side, teeth cavities face-up, onto the filter rack, into the cleaner.



### Step 2:


- Add tap water into cleaner, submerging appliance(s), below the max water level.
- Add one (1) cleaning tablet.
- Close the lid.

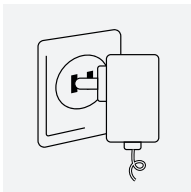


### Step 3:


- Plug in the power adapter.
- The indication light turns from black to white.

**Caution:** Do NOT run the cleaner without water added, or the cleaner will be damaged.

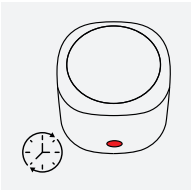

### Step 4:

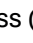- Press (  ) button. The light turns from white to orange.
- Enjoy 10 minutes of heating ( $60\pm 5^{\circ}\text{C}$ ,  $140\pm 41^{\circ}\text{F}$ ) and ultrasonic.

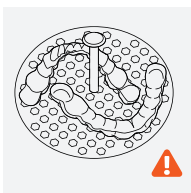


### Step 5:

- The light turns from orange to white automatically after 10 minutes.
- Unplug the cleaner and the light turns from white to black.
- Open the lid to cool down for 5 minutes.
- Lift plastic handle to take rack out.



### Step 6:

- Rinse appliance(s) with **cold** running tap water.
  - Wear freshly rejuvenated appliance(s).
- TIP:** It is normal for the appliance(s) to feel tighter immediately after rejuvenation.

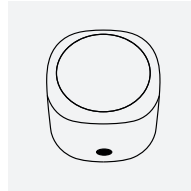

### Step 7:

- Drain the water.
- Wipe and dry the tank and the cleaner's surface.
- Store the cleaner at a cool and dry place.


**Caution:** Do NOT wash the cleaner with water to avoid electric shock and damage.

**! Caution:** Use caution and heat-resistant gloves when handling hot water. If you are under the age of 18, or have difficulty safely completing these tasks, please seek the assistance of a parent or guardian.

### TIPS

Rinse the rejuvenated appliance(s) in cold tap water and gently brush to remove any tablet residue from the surface. The residue is easily washed away with a gentle rub, leaving your appliance(s) feeling fresh and clean.

**Note:** Do NOT brush under HOT tap water, which might cause damage to the appliance(s).

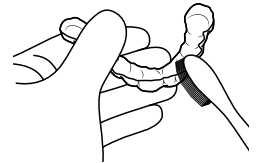

### Additional Guidance

- Make sure that your appliance(s) aren't touching each other or touching the container wall.
- HOT Surface! Be extra careful if you take out appliance(s) right after 10 minutes of rejuvenation process.

## 04 Malfunction Resolution

If there is malfunction of the ultrasonic cleaner like unable to start the machine or unable to perform rejuvenation process that you cannot handle, please contact LuxCreo CST immediately by [support@luxcreo.com](mailto:support@luxcreo.com).

### Contact and Support

Contact our Customer Success Team via:

W: [support.luxcreo.com](http://support.luxcreo.com)

E: [support@luxcreo.com](mailto:support@luxcreo.com)

T: +1 (650) 336-0888

A: 350 W. Ontario Street, Suite 700, Chicago, IL 60654